


If the connection between Epic and the ORALyzer goes down, a manual version of the OralFitnessCheck® screening will need to be ordered. After processing the sample in the ORALyzer, the result will need to be manually entered into the patient's chart.


## Order an OralFitnessCheck Screening in Epic

1. In the patient's chart, select the  Order Preference List.



2. In the search field at the top of the window, enter "oral fitness man".
3. Select Oral Fitness aMMP-8 Manually Resulted.

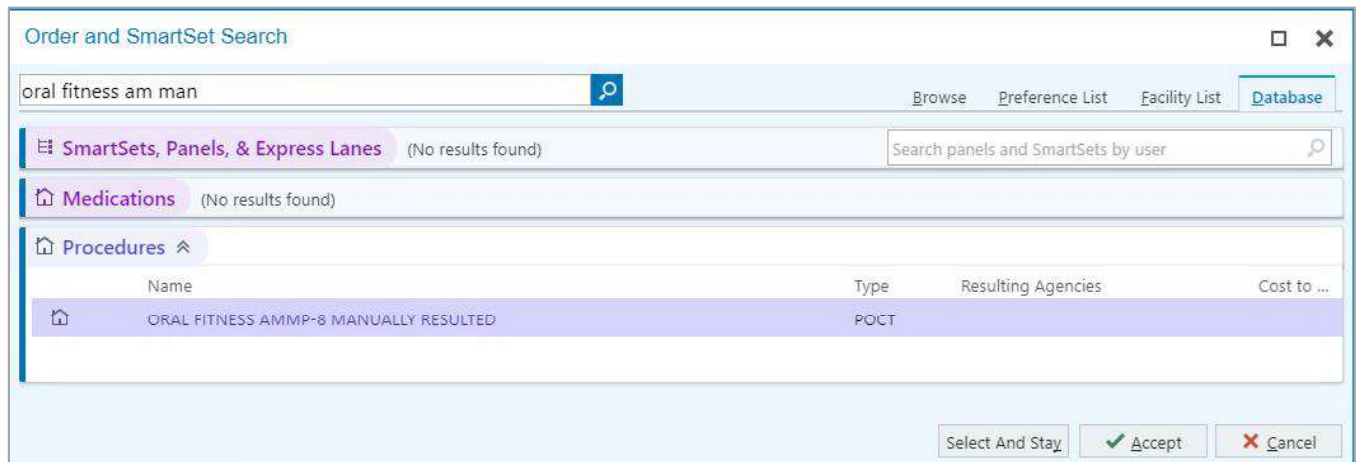


4. Click  Accept.
5. Click  Sign Orders.



If a dental assistant places an OralFitnessCheck order, the order mode will default to *Per Protocol: Cosign Required*. This order mode allows a dental assistant to sign the order. When a dental assistant signs the order, the order will transmit to the ORALyzer and the visit's dentist will receive an In Basket message to cosign the order.


## Enter OralFitnessCheck Screening Results in Epic

1. In the Treatment Plan tab, locate the POCT to Result section and click Enter Results next to the Oral Fitness AMMP-8 order.

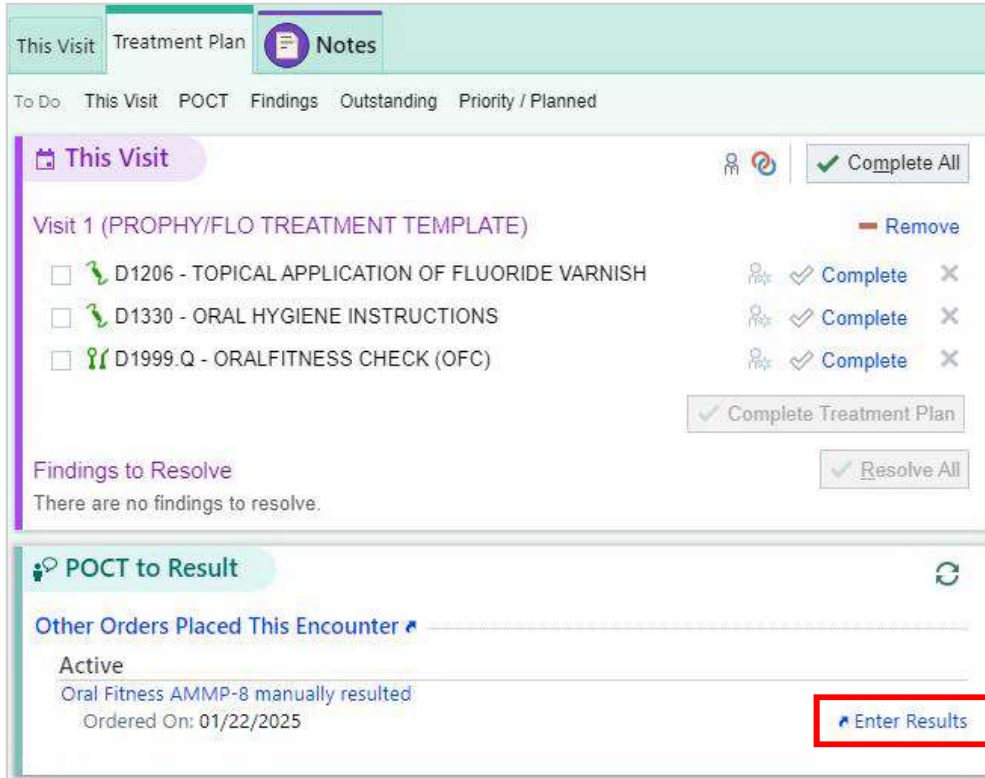

The screenshot shows the Epic Treatment Plan tab. The 'POCT to Result' section is visible, and the 'Enter Results' button is highlighted with a red box.

2. Enter the screening results, the result date, and time, and ensure the status is set to final.

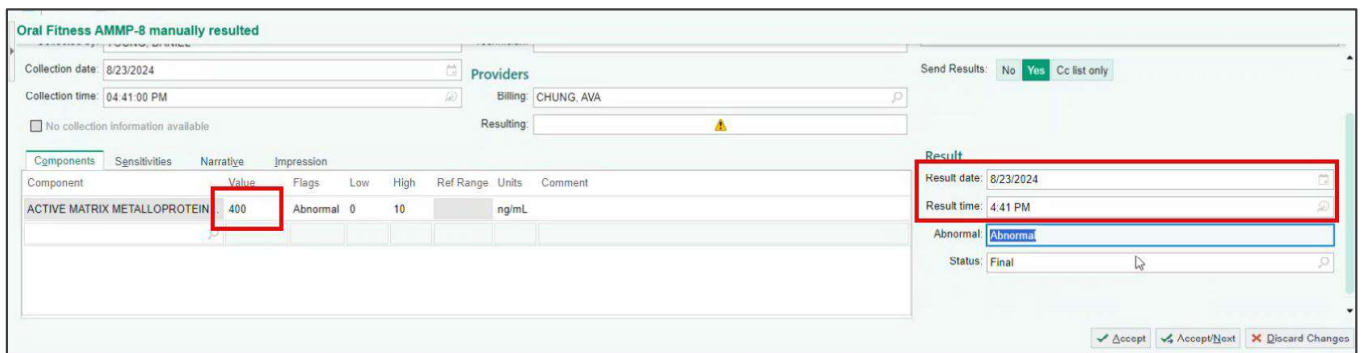

The screenshot shows the Epic Oral Fitness AMMP-8 manually resulted form. The 'Value' field (400) and the 'Result' section (Result date: 8/23/2024, Result time: 4:41 PM, Status: Final) are highlighted with red boxes.

3. Click Accept or Accept/Next.